

Histone H3 Mouse Monoclonal Antibody(Zebrafish Specific)

Catalog	TDY148C	TDY148F	Tel: 010-82908854
Quantity	50 μ L	100 μ L	Free: 400-0620-621
			Web: www.tdybio.com

For research use only.

Applications	Species Cross-Reactivity	Molecular Weight	Isotype
WB	Zebrafish	15KD	IgG1

Storage Buffer & Condition: PBS, pH 7.4, containing 0.02% **sodium azide** as Preservative and 50% Glycerol.
Store at **-20°C**. **Do not aliquot the antibody.**


Recommended dilutions: WB: 1:1,000-3,000

Optimal dilutions should be determined by the end user.

Specificity: The Histone H3 antibody can detects Zebrafish endogenous Histone H3 protein.

Alternative Names: H3 histone antibody, HIST1H3A antibody, Histone cluster 1, H3a antibody

Background: Histone H3 is one of the five main histone proteins involved in the structure of chromatin in eukaryotic cells. Core component of nucleosome. Nucleosomes wrap and compact DNA into chromatin, limiting DNA accessibility to the cellular machineries which require DNA as a template. Histones thereby play a central role in transcription regulation, DNA repair, DNA replication and chromosomal stability.

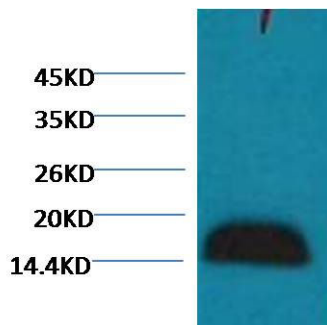

Western blot analysis of Zebrafish skeletal muscle with TDY148 diluted at 1:2,000.