

## Flotillin-1 Mouse Monoclonal Antibody(6C10)

Catalog	TDY197C	TDY197F
Quantity	50 $\mu$ L	100 $\mu$ L

Tel: 010-82908854  
Free: 400-0620-621  
Web: www.tdybio.com

**For research use only.**

Applications	Species	Cross-Reactivity	Molecular Weight	Isotype
WB		R, M	49KD	IgG1

**Storage Buffer & Condition:** PBS, pH 7.4, containing 0.02% **sodium azide** as Preservative and 50% Glycerol.  
Store at **-20°C**. **Do not aliquot the antibody.**


**Recommended dilutions: WB: 1:1,000-2,000**

**Optimal dilutions should be determined by the end user.**

**Specificity:** Antibody can detects endogenous Flotillin-1 protein.

**Alternative Names:** FLOT 1 antibody, Integral membrane component of caveolae antibody, Reggie 2 antibody

**Background:** Two flotillin members have been described, flotillin-1 and flotillin-2. May act as a scaffolding protein within caveolar membranes, functionally participating in formation of caveolae or caveolae-like vesicles.

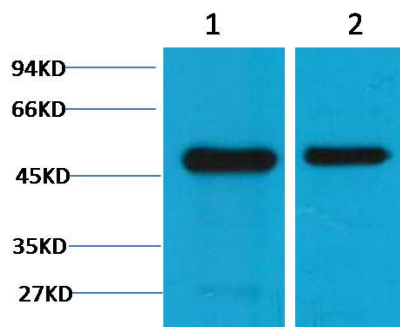

Western blot analysis of 1) Mouse Brain Tissue , 2) Rat Brain tissue with TDY197 diluted at 1:2,000.