

## Rabbit IgG(Light Chain Specific) Mouse Monoclonal Antibody(8D6)

Catalog	TDY304C	TDY304F	Tel: 010-82908854
Quantity	50 $\mu$ L	100 $\mu$ L	Free: 400-0620-621
			Web: www.tdybio.com

**For research use only.**


Applications	Species Cross-Reactivity	Molecular Weight	Isotype
WB, IHC	Rabbit	~25KD	IgG1

**Storage Buffer & Condition:** PBS, pH 7.4, containing 0.02% **sodium azide** as Preservative and 50% Glycerol.  
Store at **-20°C**. **Do not aliquot the antibody.**

**Recommended dilutions:** WB: 1:5,000-20,000 IHC: 1:200-500

**Optimal dilutions should be determined by the end user.**

**Specificity:** The antibody reacts with the light chains on rabbit IgG and with those common to other rabbit immunoglobulins. The antibody **does not react with the heavy chain of rabbit IgG**. The antibody minimal cross-reaction with **bovine, goat, Armenian hamster, horse, human, mouse, rat, and sheep** immunoglobulins, but it may cross-react with IgG from other species.

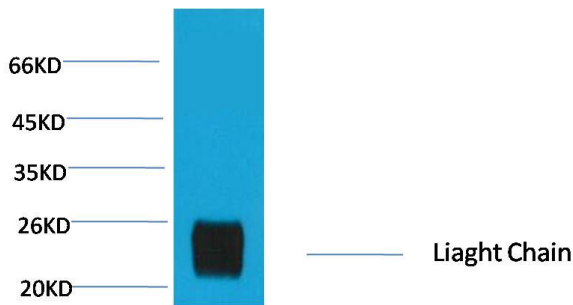

Western blot analysis of Reduced Rabbit IgG, H+L with TDY304 diluted at 1:10,000.