

Histone H3(Phospho Thr118) Rabbit Polyclonal Antibody

Catalog TDY235C TDY235F

Tel: 010-82908854

Quantity 50 μ L 100 μ L

Free: 400-0620-621

Web: www.tdybio.com

For research use only.

Applications	Species	Cross-Reactivity	Molecular Weight	Isotype
WB, IHC		H, R, M	~15-17KD	IgG

Storage Buffer & Condition: PBS, pH 7.4, containing 0.02% **sodium azide** as Preservative and 50% Glycerol.
Store at **-20°C**. **Do not aliquot the antibody.**



Recommended dilutions: WB: 1:1,000 IHC: 1:100-200

Optimal dilutions should be determined by the end user.

Specificity: The Histone H3(Phospho Thr118) Rabbit Polyclonal Antibody can detect endogenous Histone H3(Phospho Thr118) protein.

Alternative Names: H3 histone antibody, HIST1H3A antibody, Histone cluster 1, H3a antibody

Background: Histone H3 is one of the five main histone proteins involved in the structure of chromatin in eukaryotic cells. Core component of nucleosome. Nucleosomes wrap and compact DNA into chromatin, limiting DNA accessibility to the cellular machineries which require DNA as a template. Histones thereby play a central role in transcription regulation, DNA repair, DNA replication and chromosomal stability.

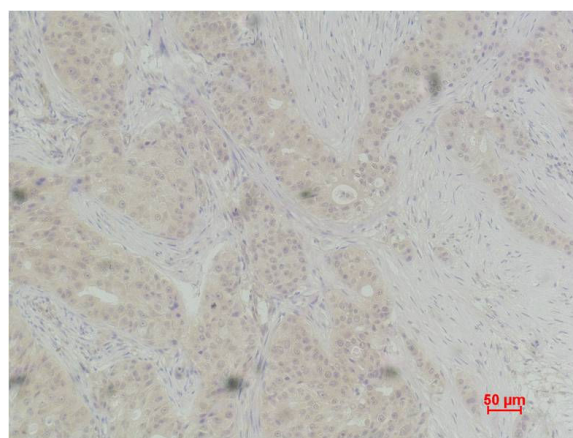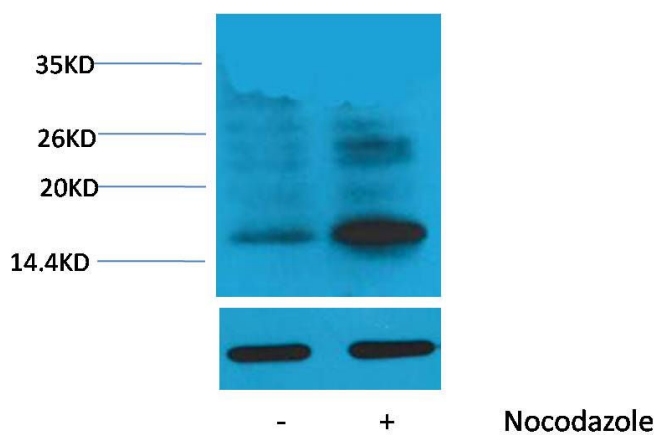

Western blot analysis of extracts from Hela cells, untreated (-) or treated with Nocodazole (1 μ g/ml, 24 hr; +), using Phospho-Histone H3 (Thr118) Rabbit pAb (upper) TDY2351:1000 and Histone H3 (1G1) Mouse mAb (lower) TDY055 1:2000.

Immunohistochemical analysis of paraffin-embedded Human Breast Carcinoma using Histone H3(Phospho Thr118) (TDY235) Rabbit pAb diluted at 1:200.