

KVβ1 Rabbit Polyclonal Antibody(A255)

Catalog TDY547C TDY547F

Tel: 010-80117836

Web: www.tdybio.com

Quantity 50μL 100μL

Entrez-Gene ID# 7811, Swiss-Prot Acc.#Q14722

For research use only.

Applications	Species Cross-Reactivity	Molecular Weight	Isotype
IHC	H,R,M	~50KD	IgG

Storage Buffer & Condition: Antigen Affinity Purified IgG in PBS, pH 7.4, containing 0.02% **sodium azide** as Preservative and 50% Glycerol.


Store at **-20°C. Do not aliquot the antibody.**

Recommended dilutions: IHC: 1:100-200


Optimal dilutions should be determined by the end user.

Specificity: Antibody can detects endogenous KVβ1 protein.

Alternative Names: hKvb3,KCNA1B,KV beta 1,Voltage-gated potassium channel beta 1 subunit

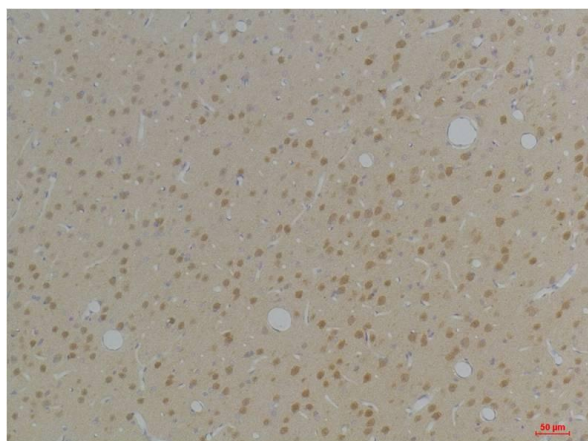

Immunohistochemical analysis of paraffin-embedded Rat BrainTissue using Kvb1 (TDY547) Rabbit pAb diluted at 1:200.

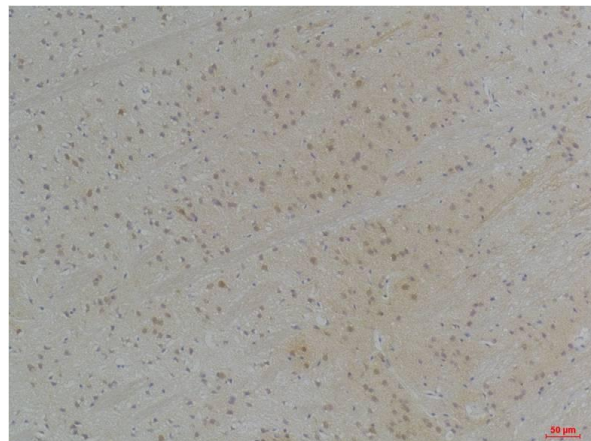

Immunohistochemical analysis of paraffin-embedded Mouse BrainTissue using Kvb1 (TDY547) Rabbit pAb diluted at 1:200.