

KV11.1 Rabbit Polyclonal Antibody(A262)

Catalog TDY553C TDY553F

Tel: 010-80117836

Web: www.tdybio.com

Quantity 50μL 100μL

Entrez-Gene ID# 169522, Swiss-Prot Acc.#Q8TDN2

For research use only.

Applications	Species Cross-Reactivity	Molecular Weight	Isotype
WB,IHC	R,M	~62KD	IgG

Storage Buffer & Condition: Antigen Affinity Purified IgG in PBS, pH 7.4, containing 0.02% **sodium azide** as Preservative and 50% Glycerol.


Store at **-20°C. Do not aliquot the antibody.**

Recommended dilutions: WB: 1:1,000-2,000 IHC: 1:100-200


Optimal dilutions should be determined by the end user.

Specificity: Antibody can detects endogenous KV11.1 protein.

Alternative Names: LQT2,HERG

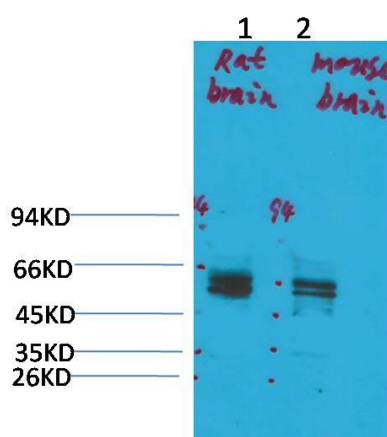

Western blot analysis of 1) Rat Brain Tissue, 2) Mouse Brain Tissue with KV11.1 Rabbit pAb TDY553 diluted at 1:2,000.

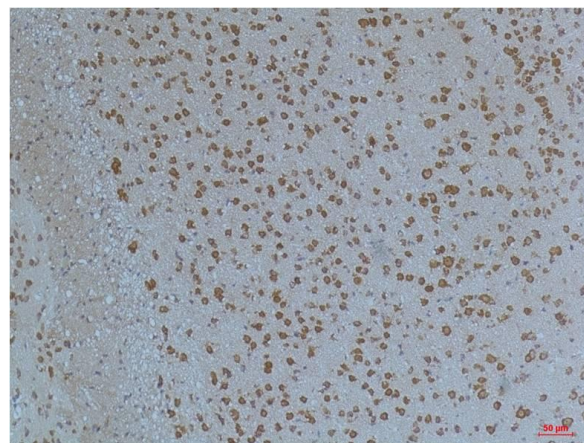

Immunohistochemical analysis of paraffin-embedded Mouse Brain Tissue using KV11.1 (TDY553) Rabbit pAb diluted at 1:200.