

KV11.3 Rabbit Polyclonal Antibody(A263)

Catalog TDY554C TDY554F

Tel: 010-80117836

Web: www.tdybio.com

Quantity 50μL 100μL

Entrez-Gene ID# 90134, Swiss-Prot Acc.#Q9NS40

For research use only.

Applications	Species Cross-Reactivity	Molecular Weight	Isotype
WB,IHC	H,R,M	~83,135KD	IgG

Storage Buffer & Condition: Antigen Affinity Purified IgG in PBS, pH 7.4, containing 0.02% **sodium azide** as Preservative and 50% Glycerol.


Store at **-20°C. Do not aliquot the antibody.**

Recommended dilutions: WB: 1:1,000-2,000 IHC: 1:100-200


Optimal dilutions should be determined by the end user.

Specificity: Antibody can detects endogenous KV11.3 protein.


Alternative Names: EAG2,H EAG2,KCNH5

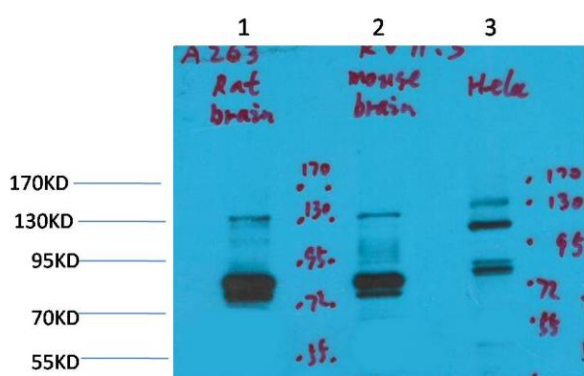

Western blot analysis of 1) Rat Brain Tissue, 2) Mouse Brain Tissue, 3) HeLa with KV11.3 Rabbit pAb TDY554 diluted at 1:2,000.

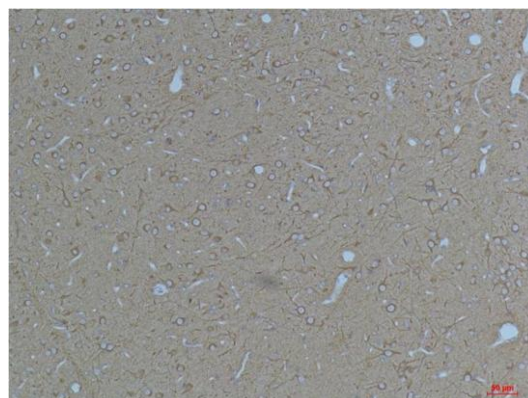

Immunohistochemical analysis of paraffin-embedded Rat Brain Tissue using KV11.3 (TDY554) Rabbit pAb diluted at 1:200.

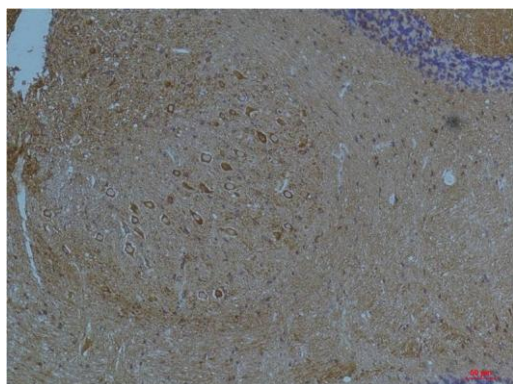

Immunohistochemical analysis of paraffin-embedded Mouse Brain Tissue using KV11.3 (TDY554) Rabbit pAb diluted at 1:200.

Applications: WB-Western blot IHC-Immunochemistry IF-Immunofluorescence IP-Immunoprecipitation ChIP-Chormatin Immunoprecipitation
Reactivity: H-Human R-Rat M-Mouse Mk-Monkey Dg-Dog Ch-Chicken Hm-Hamster Rb-Rabbit Sh-Sheep Pg-Pig Z-Zebrafish