

PPAR Delta Mouse Monoclonal Antibody(2F9)

Catalog TDY569C TDY569F

Tel: 010-80117836

Web: www.tdybio.com

Quantity 50μL 100μL

Entrez-Gene ID#5469 , Swiss-Prot Acc.#Q03181

For research use only.

Applications	Species Cross-Reactivity	Molecular Weight	Isotype
WB,IHC	H,R,M	~50KD	IgG1

Storage Buffer & Condition: PBS, pH 7.4, containing 0.02% **sodium azide** as Preservative and 50% Glycerol.

Store at **-20°C. Do not aliquot the antibody.**


Recommended dilutions: WB: 1:1,000-2,000 IHC: 1:100-200

Optimal dilutions should be determined by the end user.


Specificity: The PPAR Delta antibody can detects endogenous PPAR Delta protein.

Alternative Names: FAAR, NR1C2, NUC1, Peroxisome proliferative activated receptor delta

Background:

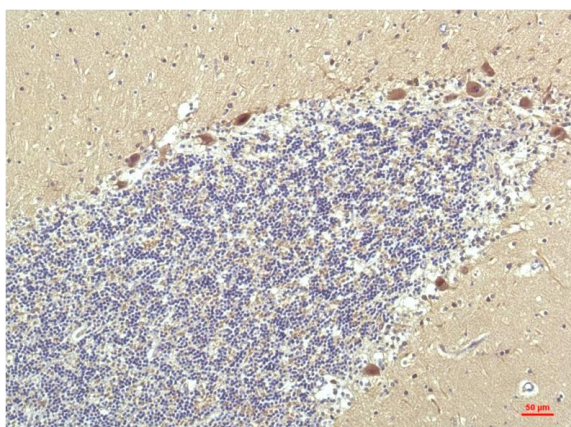

Immunohistochemical analysis of paraffin-embedded Human Brain Tissue using PPAR Delta (TDY569) Mouse mAb diluted at 1:200.

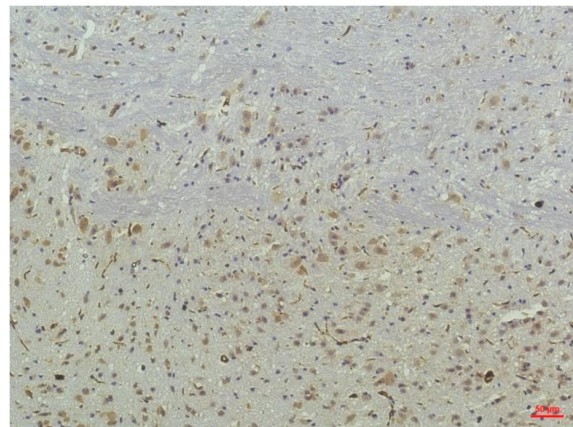

Immunohistochemical analysis of paraffin-embedded Mouse Brain Tissue using PPAR Delta (TDY569) Mouse mAb diluted at 1:200.