

Tau Mouse Monoclonal Antibody(10E3)

Catalog TDY640C TDY640F

Tel: 010-80117836

Web: www.tdybio.com

Quantity 50μL 100μL

Entrez-Gene ID:4137 , Swiss-Prot Acc.P10636

For research use only.

Applications	Species Cross-Reactivity	Molecular Weight	Isotype
WB,IHC	H,R,M	~50-80KD	IgG1

Storage Buffer & Condition: PBS, pH 7.4, containing 0.02% **sodium azide** as Preservative and 50% Glycerol.

Store at **-20°C. Do not aliquot the antibody.**


Recommended dilutions: WB: 1:1,000-2,000 IHC: 1:100-200

Optimal dilutions should be determined by the end user.


Specificity: Antibody can detects endogenous Tau protein.

Alternative Names: DDPAC,MAPTL,MGC134287,PHF,PPND

Background: Tau is a heterogeneous microtubule-associated protein that promotes and stabilizes microtubule assembly, especially in axons.

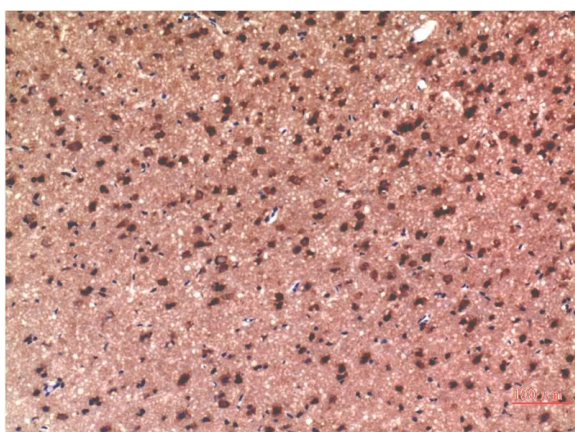

Immunohistochemical analysis of paraffin-embedded Rat Brain Tissue using Tau (TDY640) Mouse mAb diluted at 1:200.

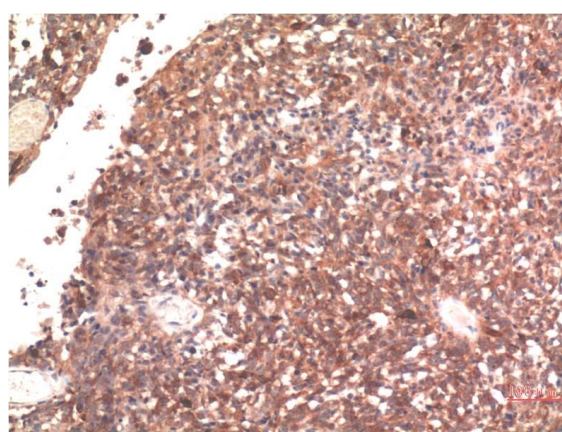

Immunohistochemical analysis of paraffin-embedded Human Brain Tissue using Tau (TDY640) Mouse mAb diluted at 1:200.