

VEGF Rabbit Polyclonal Antibody(VEGF)

Catalog TDY649C TDY649F

Tel: 010-80117836

Web: www.tdybio.com

Quantity 50μL 100μL

Entrez-Gene ID:7422 , Swiss-Prot Acc.P15692

For research use only.

Applications	Species Cross-Reactivity	Molecular Weight	Isotype
IHC	H,R,M	~23,45KD	IgG

Storage Buffer & Condition: Antigen Affinity Purified IgG in PBS, pH 7.4, containing 0.02% **sodium azide** as Preservative and 50% Glycerol.


Store at **-20°C**. **Do not aliquot the antibody.**

Recommended dilutions: IHC: 1:200-500


Optimal dilutions should be determined by the end user.

Specificity: Antibody can detects endogenous VEGF protein.


Alternative Names: VEGFA antibody, VPF antibody, VEGF120 antibody, Vascular endothelial growth factor antibody

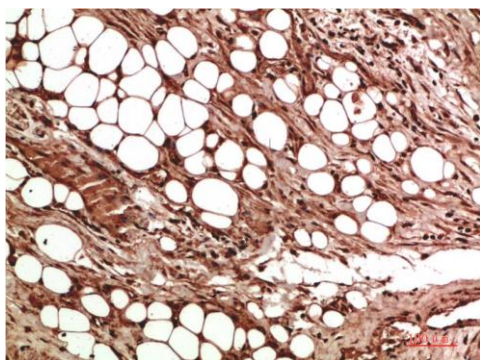

Immunohistochemical analysis of paraffin-embedded Human Liver Carcinoma Tissue using VEGF (TDY649) Rabbit pAb diluted at 1:500.

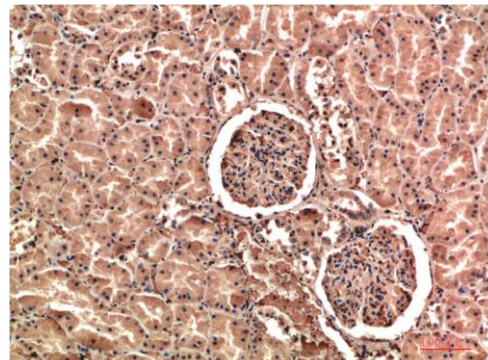

Immunohistochemical analysis of paraffin-embedded Human Kidney Tissue using VEGF (TDY649) Rabbit pAb diluted at 1:500.

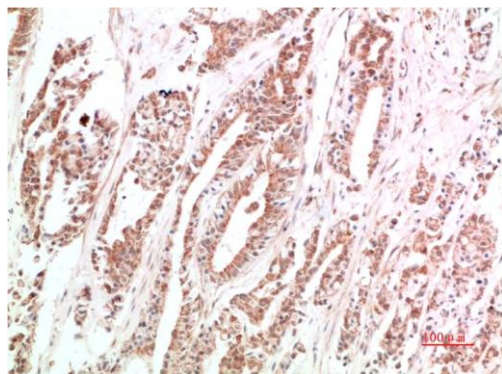

Immunohistochemical analysis of paraffin-embedded Human Stomach Carcinoma Tissue using VEGF (TDY649) Rabbit pAb diluted at 1:500.