

JAK2 Rabbit Polyclonal Antibody(F175)

Catalog TDY688C TDY688F

Tel: 010-80117836

Web: www.tdybio.com

Quantity 50μL 100μL

Entrez-Gene ID:3717 , Swiss-Prot Acc.O60674

For research use only.

Applications	Species Cross-Reactivity	Molecular Weight	Isotype
IHC	H,R,M	~120KD	IgG

Storage Buffer & Condition: Antigen Affinity Purified IgG in PBS, pH 7.4, containing 0.02% **sodium azide** as Preservative and 50% Glycerol.

Store at **-20°C**. **Do not aliquot the antibody.**


Recommended dilutions: IHC: 1:100-200

Optimal dilutions should be determined by the end user.


Specificity: The Antibody can detects endogenous JAK2 proteins.

Alternative Names: Janus kinase 2,JTK10,Tyrosine protein kinase Jak2, THCYT3

Background: Members of the Janus family of tyrosine kinases (Jak1, Jak2, Jak3 and Tyk2) are activated by ligands binding to a number of associated cytokine receptors.

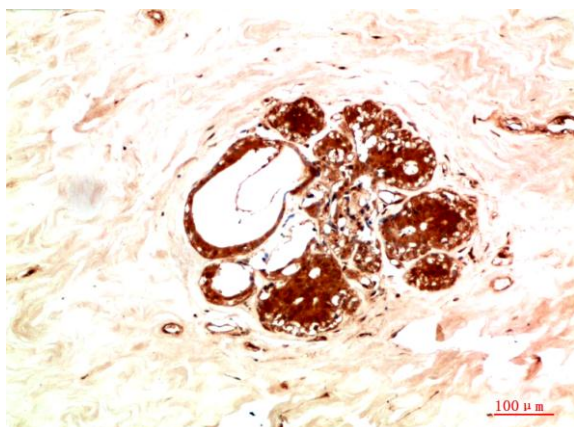

Immunohistochemical analysis of paraffin-embedded Human Breast Carcinoma Tissue using JAK2 (TDY688) Rabbit pAb diluted at 1:200

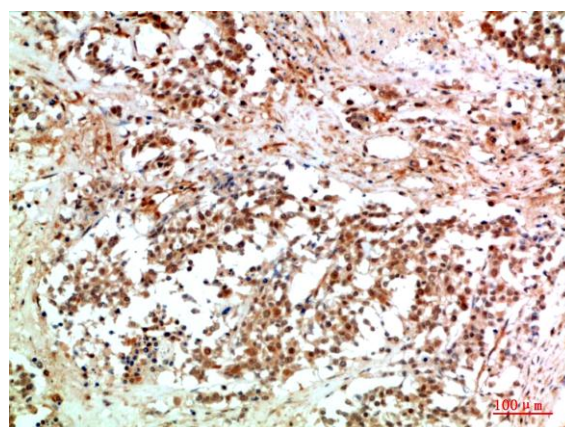

Immunohistochemical analysis of paraffin-embedded Human Lung Carcinoma Tissue using JAK2 (TDY688) Rabbit pAb diluted at 1:200