

HIF-1 β /ARNT Mouse Monoclonal Antibody(4D9)

Catalog TDY1096C TDY1096F

Tel: 010-80117836

Web: www.tdybio.com

Quantity 50 μ L 100 μ L

Entrez-Gene ID#405 , Swiss-Prot Acc.#P27540

For research use only.

Applications	Species Cross-Reactivity	Molecular Weight	Isotype
WB	H, M, R	~87KD	IgG1

Storage Buffer & Condition: PBS, pH 7.4, containing 0.02% **sodium azide** as Preservative and 50% Glycerol.
Store at **-20°C. Do not aliquot the antibody.**


Recommended dilutions: WB: 1:1,000-2,000

Optimal dilutions should be determined by the end user.

Specificity: Antibody can detects endogenous Human HIF-1 β /ARNT protein.

Alternative Names: HIF 1 beta, ARNT, Tango, hypoxia inducible factor 1 beta, Nrnt, Dioxin receptor, bHLH2e

Background: Hypoxia-inducible factors (HIFs) are transcription factors that respond to changes in available oxygen in the cellular environment, to be specific, to *decreases* in oxygen, or hypoxia. This protein is required for the ligand-binding subunit to translocate from the cytosol to the nucleus after ligand binding.

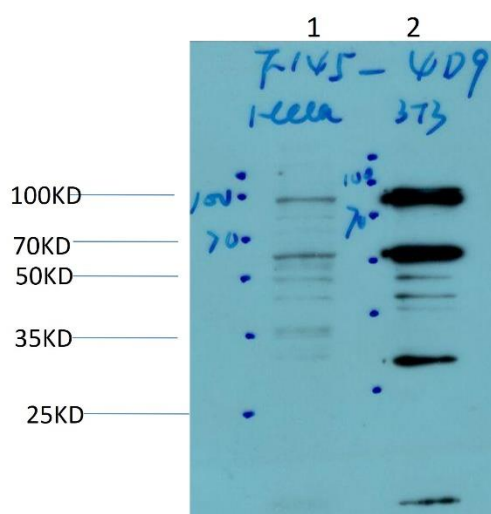

Western blot analysis of 1)HeLa Cell, 2) 3T3 Cell Lysate using HIF-1 β /ARNT

(TDY1096) Mouse Monoclonal Antibody diluted at 1:2,000.