

c-Jun Mouse Monoclonal Antibody(4H7)

Catalog	TDY1092C	TDY1092F	Tel: 010-82908854
			Web: www.tdybio.com
Quantity	50μL	100μL	Entrez-Gene ID#3725, Swiss-Prot Acc.#P05427

For research use only.

Applications	Species Cross-Reactivity	Molecular Weight	Isotype
WB	M	39-43KD	IgG1

Storage Buffer & Condition: PBS, pH 7.4, containing 0.02% **sodium azide** as Preservative and 50% Glycerol.

Store at **-20°C. Do not aliquot the antibody.**


Recommended dilutions: WB: 1:2,000-5,000

Optimal dilutions should be determined by the end user.

Specificity: Antibody can detects endogenous c-Jun proteins.

Alternative Names: Activator Protein 1 antibody, AP1 antibody, p39 antibody

Background: c-Jun is a member of the Jun Family containing c-Jun, JunB and JunD, and is a component of the transcription factor AP-1 (activator protein-1).

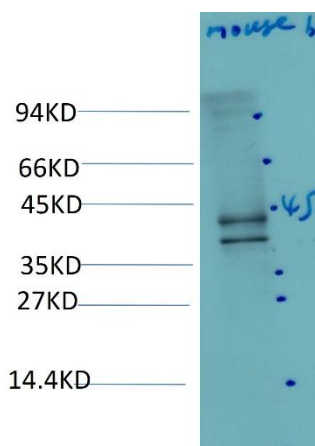

Western blot analysis of Mouse Brain Tissue Lysate using c-Jun (TDY1092)

Mouse Monoclonal Ab diluted at 1:5,000.